

ALCONOX PRODUCT LIST

Available worldwide from leading laboratories, scientific, healthcare and industrial suppliers.


Biodegradable. Concentrated. Interfering residue free.

Alconox Inc. detergents set the standard for critical cleaning. Scientists, engineers and professionals specify Alconox Inc. detergents in their cleaning procedures for optimum economy, environmental safety and performance.







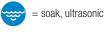





ALCONOX®

Powdered Precision Cleaner Original general-purpose powder for manual, soak, and ultrasonic cleaning of glass, metal, and plastic. No interfering residues. Has inhibitory residue test. USDA/

FDA compliant. Dilute 1:100. pH 9.5

| Size | Mfr. Cat No | |
|-------------------------|-------------|--|
| 4 lb box (1.8 kg) | 1104-1 | |
| CS 9 x 4 lb (9 x 1.8 kg |) 1104 | |
| 25 lb box (11 kg) | 1125 | |
| 50 lb box (22.5 kg) | 1150 | |
| 100 lb drum (45 kg) | 1101 | |
| 300 lb drum (136 kg) | 1103 | |
| 50 x .5 oz packs (50 x | 15g) 1112-1 | |
| CS 12 x 50 packs | 1112 | |



LIQUINOX® Critical Cleaning Liquid Detergent

Original general-purpose, phosphate-free liquid for manual, soak, and ultrasonic cleaning of glass, metal, and plastic. No

interfering residues. Has inhibitory residue test. USDA/ FDA compliant, Dilute 1:100, pH 8.5

| 1 qt bottle (.95 L) | 1232-1 | |
|---------------------------|--------|--|
| CS 12 x 1 qt (12 x .95 L) | 1232 | |
| 1 gal bottle (3.8 L) | 1201-1 | |
| CS 4 x 1 gal (4 x 3.8 L) | 1201 | |
| 5 gal jerrycan (19 L) | 1205 | |
| 15 gal drum (57 L) | 1215 | |
| 55 gal drum (209 L) | 1255 | |



CITRANOX®

Liquid Acid Cleaner and Detergent

Liquid acid cleaner for manual, soak and ultrasonic cleaning or acid neutralizing. Corrosion-inhibited for glassware, metals,

and plastic. Surfactant rinse aids for no interfering residues. Has inhibitory residue test. USDA and FDA compliant. Dilute 1:100. pH 2.5

| Size | Mfr. Cat No |
|--------------------------|-------------|
| 1 qt bottle (.95 L) | 1832-1 |
| CS 12 x 1 qt (12 x .95 | L) 1832 |
| 1 gal bottle (3.8 L) | 1801-1 |
| CS 4 x 1 gal (4 x 3.8 L) | 1801 |
| 5 gal jerrycan (19 L) | 1805 |
| 15 gal drum (57 L) | 1815 |
| 55 gal drum (209 L) | 1855 |
| | |



DETONOX®

Ultimate Precision Cleaner Non-caustic, liquid concentrate

detergent for hand and ultrasonic use on challenging residues. Exceptional removal of cosmetic and pharmaceutical polymers,

sticky extracts, resins, creams and lotions. Concentrated, highly emulsifying and penetrating formula is free rinsing for reliable results, and leaves no interfering residues. Dilute 1:100, pH 10.5

| 1 gal bottle | 2301-1 | |
|--------------------------|--------|--|
| CS 4 x 1 gal (4 x 3.8 L) | 2301 | |
| 5 gal jerrycan (19 L) | 2305 | |
| 55 gal drum (208 L) | 2355 | |
| | | |



(1)

LUMINOX®

Low-Foaming Neutral pH Liquid Detergent

Low-foaming, neutral pH detergent for use in manual, ultrasonic, and machine cleaning. Safe-handling detergent used

for noncorrosive, light-duty cleaning. No chelation or alkaline waste treat problems. Free rinsing, eliminating interfering residues. Dilute 3:100. pH 7.0

| Size I | Mfr. Cat No | |
|--------------------------|-------------|--|
| 1 gal bottle (3.8 L) | 1901-1 | |
| CS 4 x 1 gal (4 x 3.8 L) | 1901 | |
| 5 gal jerrycan (19 L) | 1905 | |
| 15 gal drum (57 L) | 1915 | |
| 55 gal drum (209 L) | 1955 | |
| | | |



TERGAZYME®

Enzyme-Active Powdered Detergent

Enzymatic cleaner for manual, soak, and ultrasonic cleaning of proteins, tissue, and body fluids from glass, metals,

and plastic, with no interfering residues. Cleans filter membranes, fermenters and bioreactors. Has inhibitory residue test. FDA/USDA compliant. Dilute 1:100. pH 9.5


| 4 lb box (1.8 kg) | 1304-1 | |
|--------------------------|--------|--|
| CS 9 x 4 lb (9 x 1.8 kg) | 1304 | |
| 25 lb box (11 kg) | 1325 | |
| 50 lb box (22.5 kg) | 1350 | |
| 100 lb drum (45 kg) | 1301 | |
| 300 lb drum (136 kg) | 1303 | |









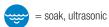



CITRAJET® **Low-Foam Liquid Acid** Cleaner/Rinse

Liquid acid cleaner for machine, spray-CIP, manual, and ultrasonic cleaning or acid neutralizing. Corrosion-inhibited

for glassware, metals, and plastic. Surfactant rinse for no interfering residues. Has inhibitory residue test. FDA compliant, Dilute 1:100, pH 2.5

| Size | Mfr. Cat No |
|--------------------------|----------------|
| 1 gal bottle (3.8 L) | 2001-1 |
| CS 4 x 1 gal (4 x 3.8 L) | 2001 |
| 5 gal can (19 L) | 2005 |
| 15 gal drum (57 L) | 2015 |
| 30 gal drum (114 L) | 2030 |
| 55 gal drum (209 L) | 2055 |



DETOJET® Low-Foaming Liquid Detergent

Heavy-duty alkaline degreaser. Nonionic liquid for labware and machine washers, manual, and

ultrasonic cleaning. Corrosion-inhibited for glass, metals, and plastic. No interfering residues. Has inhibitory residue test. USDA/FDA compliant. Dilute 1:200. pH 12

| 1qt bottle (.95 L) | 1632-1 |
|---------------------------|--------|
| CS 12 x 1 qt (12 x .95 L) | 1632 |
| 1gal bottle (3.8 L) | 1601-1 |
| CS 4 x 1 gal (4 x 3.8 L) | 1601 |
| 5 gal can (19 L) - Air OK | 1605 |
| 15 gal drum (57 L) | 1615 |
| 30 gal drum (114 L) | 1630 |
| 55 gal drum (209 L) | 1655 |
| | |



ALCOJET® Low-Foaming Powdered Detergent

General-purpose nonionic powder for labware and machine washers, manual, and ultrasonic cleaning.

Corrosion-inhibited for glassware, metals, and plastic. No interfering residues. Has inhibitory residue test. USDA authorized, Dilute 1:200, pH 11.5

| Size | Mfr. Cat No | |
|--------------------------|----------------|--|
| 4 lb box (1.8 kg) | 1404-1 | |
| CS 9 x 4 lb (9 x 1.8 kg) | 1404 | |
| 25 lb box (11 kg) | 1425 | |
| 50 lb box (22.5 kg) | 1450 | |
| 100 lb drum (45 kg) | 1401 | |
| 300 lb drum (136 kg) | 1403 | |
| | | |



ALNOCHROMIX^{TO}

Oxidizing Acid Additive for Glass Cleaning

Concentrated inorganic persulfate based additive to sulfuric acid to make

a metal-free and surfactant-free glassware cleaner to be used in tissue culture, trace metal, and enzyme kinetics research, etc. Nonmetallic, non-carcinogenic, and easier to dispose replacement for chromic acid.

| 1 box 10 x 3.1 oz (87 g) pouches | 2510 |
|--|------|
| 4 boxes of 10 x 3.1 oz (87 g) pouches | 2540 |
| 12 boxes of 10 x 3.1 oz (87 g) pouches | 2512 |
| 1 pair silver gloves | 2521 |
| 3 safety vent caps | 2513 |
| Kit 3 vent cap+3 pair gloves | 2533 |



TERGAJET®

Low-Foaming Phosphate-Free Powder

Phosphate-free, nonionic powder for machine washers, hand, and ultrasonic cleaning. Corrosion-inhibited for

labware. For water, microbiology, and environmental labs. No interfering residues. Has inhibitory residue test. Dilute 1:200. pH 11.5

| Size | Mfr. Cat No | |
|--------------------------|----------------|--|
| 4 lb box (1.8 kg) | 2204-1 | |
| CS 9 x 4 lb (9 x 1.8 kg) | 2204 | |
| 25 lb box (11 kg) | 2225 | |
| 50 lb box (22.5 kg) | 2250 | |
| 100 lb drum (45 kg) | 2201 | |
| 300 lb drum (136 kg) | 2203 | |
| | | |



SOLUJET®

Low-Foaming Phosphate-Free Liquid

Nonionic liquid for machine, spray-CIP, manual, soak, and ultrasonic cleaning. High emulsifying and penetrating

formula. Corrosion-inhibited for glass, metals, and plastic. No interfering residues. Has inhibitory residue test, FDA compliant, Dilute 1:200, pH 12

| | F |
|---------------------------|--------|
| 1 gal bottle (3.8 L) | 2101-1 |
| CS 4 x 1 gal (4 x 3.8 L) | 2101 |
| 5 gal can (19 L) - Air OK | 2105 |
| 15 gal drum (57 L) | 2115 |
| 30 gal drum (114 L) | 2130 |
| 55 gal drum (209 L) | 2155 |



DETERGENT 8° **Low-Foaming Ion-Free Liquid Detergent**

Nonionic detergent for PC board washing machines, manual, and ultrasonic cleaning. Ion- and chelate-

free for no conductive residues on electronics and no chelate interference in waste. Phosphate-free. Dilute 3:100. pH 11

| Size | Mfr. Cat No |
|--------------------------------|----------------|
| 1 gal bottle (3.8 L) | 1701-1 |
| CS 4 x 1 gal (4 x 3.8 L) | 1701 |
| 5 gal jerrycan (19 L) - Air Ok | (1705 |
| 15 gal drum (57 L) | 1715 |
| 55 gal drum (209 L) | 1755 |



KEYLAJET®

Low-Foaming High Alkaline Liquid

Corrosion inhibited formulation recommended for hard metals, steels.

stainless steel, porcelain, ceramic, rigid plastics, cement, and elastomers. Free rinsing to give you reliable results and no interfering residues. Dilute 1:200. pH 13+

| 5 gal can (19 L) - Air OK | 2405 |
|---------------------------|-------|
| 15 gal drum (57 L) | 2415 |
| 55 gal drum (208 L) | 2455 |
| 265 gal tote (1003 L) | 24265 |



ALCOTABS®

Tablet Pipet Detergent Effervescent detergent tablet for maximum cleaning in siphon pipet washers. Slow release bubbles and free rinsing give no interfering

residues. Has inhibitory residue test. Use 2 tabs/wash.

| 100 tablet bottle | 1500-1 |
|--------------------|--------|
| CS 6 x 100 bottles | 1500 |

